


Autosampler Syringes



Ordering Information				
Volume, Code, Gauge, Length, Point Style	QTY	OEM No.	Thermo No.	Part No.
Hamilton Syringes				
10µL, N, 26s, 80mm, Cone for AS 200, TriPlus	1	202066	36502019	24-202066
5µL, N, 26s, 50mm, Cone for AS 2000/3000	1	204051	36500505	24-204051
10µL, N, 26s, 50mm, Cone for AS 2000/3000	1	204052	365D3711	24-204052
SGE Syringes				
10µL, F, 26s, 80mm, Cone for AS 200, TriPlus	1	002992	36502019	21-002992
10µL, R, 26s, 80mm, Cone for AS 200/800	1	002993		21-002993
10µL, F, 26s, 50mm, Cone for AS 2000/3000	1	002980	365D3711	21-002980
10µL, R, 26s, 50mm, Cone for AS 2000/3000	1	002982	365D1841	21-002982
Needles 10µL, 26s, 80mm, Cone	1	031535		21-031535
Needles 10µL, 26s, 50mm, Cone	1	037785		21-037785

Needle Termination Codes	
Hamilton	
N	Cemented Needle
SGE	
F	Fixed Needle
R	Removable Needle

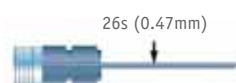

Point Style Cone
Cone tip

Septum Nut



Septum nut provides consistent septum tightness, fewer bent needles, and easy leak detection. The nut contains a needle guide, ensuring that the needle consistently penetrates the septum in the same place, prolonging septum life. The guide also prevents the needle from striking the edge of the column or bending during insertion.

Septum Nuts for Thermo Trace GCs

- Ensure a leak-tight injection port, increase septum lifetime, and decrease maintenance requirements.
- Needle guide allows easy penetration and prevents septum coring.
- Manual injection septum nut allows use of 23-gauge needles for on-column injections.
- Made of high-quality stainless steel.
- Include septum port and holder.

Ordering Information			
Description	QTY	Thermo No.	Part No.
Septum Nut for Varian			
Septum Nut (for 23-Gauge Needles)	1		065B-0011